


Agilon™

COLLAGEN ENHANCED

Based on our flagship product, Morpheus™, **Agilon moldable collagen enhanced bone graft substitute** offers all of the same handling and healing benefits, now with the addition of bovine collagen. Agilon moldable consists of 1-2mm proprietary TrelCor™ granules combined with a rapidly absorbing, biocompatible organic binder and blended with type 1 collagen. The added collagen fibers increase overall cohesiveness, providing a firmer feel without sacrificing moldability. Additionally, the added collagen minimizes residual material left on surgical gloves, aiding in optimal handling characteristics during mixing and intra-operative placement. The organic binder resists breakdown during handling and placement, then is rapidly absorbed, leaving behind the porous TrelCor granules.

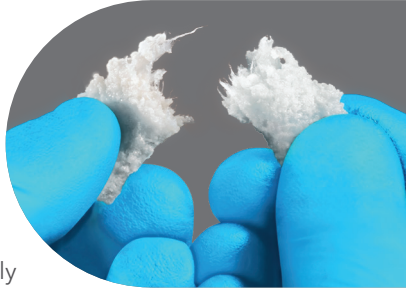

Configuration: Putty

Chemistry: Calcium phosphate/calcium carbonate (CP/CC) TrelCor granules and type 1 bovine collagen in an organic binder

Architecture: Interconnected porosity

Granule Pore Size: 500 microns

Binder Resorption Profile: 1-2 days

Storage and Handling: Room temperature

Agilon DIFFERENCE

HANDLING

- Very firm but moldable and mixable
- Stays in place and resists migration and irrigation

CHEMISTRY

- Resorbable, biocompatible, CP/CC TrelCor granules
- Organic binding agent and type 1 collagen

ARCHITECTURE

- Full interconnected porosity of TrelCor granules
- Putty-like consistency of binder, collagen and granule composite

Agilon™

GROW WITH US

At Biogenix we recognize that without great sales and distribution partners, our innovative products will not reach the surgeon or patient. This is why we prioritize and nurture our relationships with our distributors. Experience the difference a true partnership can make in your business. Contact us about distribution opportunities.

Catalog Number	Volume
----------------	--------

012-AGL	12cc
---------	------

006-AGL	6cc
---------	-----

003-AGL	3cc
---------	-----

1641 McGaw Avenue Irvine, CA, 92614
T 949-253-0994 F 949-266-5800 biogenix.com