


Agilon

MOLDABLE

COLLAGEN ENHANCED

Based on our flagship product, Morpheus[®], **Agilon Moldable is an advanced bone graft** that offers all of the same healing benefits, now with even better handling due to the addition of collagen. Developed with Biogenix proprietary TrelCor[™] technology, Agilon Moldable consists of 1-2mm TrelCor granules suspended in a rapidly absorbing, biocompatible organic binder and blended with type 1 collagen. The added collagen fibers increase overall cohesiveness, providing a firmer feel without sacrificing moldability. Additionally, the added collagen minimizes residual material left on surgical gloves, aiding in optimal handling during mixing and intra-operative placement. The organic binder resists breakdown during handling and placement, and is rapidly absorbed, leaving behind the porous TrelCor granules.

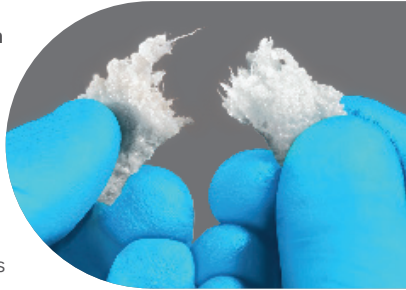

Configuration: Moldable

Chemistry: Nanocrystalline hydroxycarbanoapatite (HCA)/ calcium carbonate (CC) TrelCor granules and type 1 bovine collagen in an organic binder

Architecture: Biomimetic structure with interconnected porosity

Granule Pore Size: 500 microns

Binder Absorption Profile: 1-2 days

Granule Resorption Profile: 6-12 months

Storage and Handling: Room temperature; no reconstitution required

Agilon DIFFERENCE

TRELCOR TECHNOLOGY

- Nanocrystalline structure provides an optimal surface for cell attachment and is actively involved in bone regeneration
- HCA surface composition closely mimics bone mineral and accelerates the healing process

ARCHITECTURE

- Biomimetic interconnected pore structure is similar to cancellous bone and allows for bone growth throughout the entire graft
- High granule content ensures direct contact with surrounding bone

BEST-IN-CLASS HANDLING

- Firm consistency that is moldable, cohesive, and mixable with autograft
- Stays in place, holds its shape, and resists irrigation
- No hydration required

Agilon[®]

MOLDABLE

COLLAGEN ENHANCED
BONE GRAFT SUBSTITUTE

GROW WITH US

At Biogenix we recognize that without great sales and distribution partners, our innovative products will not reach the surgeon or patient. This is why we prioritize and nurture our relationships with our distributors. Experience the difference a true partnership can make in your business. Contact us about distribution opportunities.

Catalog Number	Volume
----------------	--------

012-AGL	12cc
---------	------

006-AGL	6cc
---------	-----

003-AGL	3cc
---------	-----

1641 McGaw Avenue Irvine, CA, 92614
T 949-253-0994 F 949-266-5800 biogenix.com

BIOGENNIX SHAPING FUSION